


USE OF THE BREEZHALER®



Drug : _____
 Indication (maintenance or rescue or maintenance+rescue): _____
 Number of inhalations: _____ Number of times per day: _____
 Other information: _____

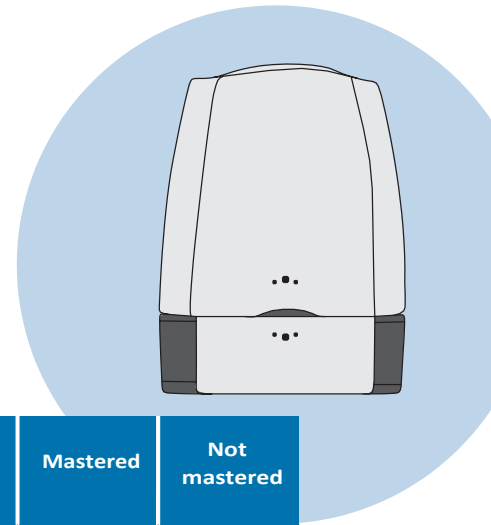

		Mastered	Not mastered
1	Prepare the device by removing the protective cap		
2	Hold the base firmly and tip the mouthpiece to open the inhaler		
3	Remove one capsule from the blister pack immediately before use: a. lift the film and release a single capsule		
4	Insert the capsule in the compartment provided for this purpose		
5	Close the mouthpiece until it clicks		
6	Pierce the capsule by holding the inhaler upwards and pressing the 2 buttons simultaneously to the bottom, once (you should hear a click), then release the buttons		
7	Breathe out normally next to the device, emptying the lungs without forcing		
8	Place the mouthpiece between the lips and teeth: a. The device is horizontal and the buttons should be on the left and right		
9	Inhale by mouth quickly and as deeply as possible: a. You should hear a whirring sound. If not, replace the capsule by tapping the base of the device		
10	Move the device away from the mouth		
11	Take a break: hold your breath for 5 to 10 seconds		
12	Breathe out through the nose		
13	Open the mouthpiece, check that there is no powder left in the capsule. If any remains, close the cap and repeat steps 8 to 13		
14	Throw the capsule away without touching it with your hands		

REMARKS

To limit the adverse effects of the medication, rinse your mouth with water (do not swallow the water)

MAINTENANCE:



RQESR

RESEAU QUÉBÉCOIS D'ÉDUCATION
 EN SANTÉ RESPIRATOIRE
 La référence des professionnels