


USE OF THE INHUB[®]



Drug : _____
 Indication (maintenance or rescue or maintenance+rescue): _____
 Number of inhalations: _____ Number of times per day: _____
 Other information: _____

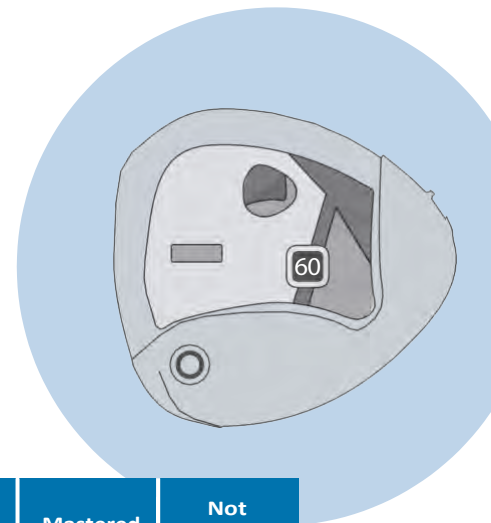

		Mastered	Not mastered
1	Remove the device from the foil package at the time of use and write the expiration date on the foil package (30 days after the opening date) : a. The device is ready to use		
2	Hold the inhaler upright in one hand and with the other hand on the notch, lower the mouthpiece cap		
3	Holding the inhaler in an upright position with one hand, push the yellow lever down with the other hand to the end, a click could be heard: a. The dose is ready to be inhaled		
4	Keep the device in an upright position		
5	Breathe out normally next to the device, emptying the lungs without forcing		
6	Place the mouthpiece in the mouth (between the lips and teeth)		
7	Breathe in through the mouth quickly and deeply: a. Be careful not to obstruct the vents		
8	Move the device away from the mouth		
9	Respect a pause: hold your breath for 5 to 10 seconds		
10	Breathe out through the nose		
11	Close the cap, the yellow lever will automatically snap back into place		

REMARKS

To limit the adverse effects of the medication, rinse your mouth with water (do not swallow the water).

MAINTENANCE:

HOW TO KNOW IF IT IS EMPTY:



RQESR

RÉSEAU QUÉBÉCOIS D'ÉDUCATION
EN SANTÉ RESPIRATOIRE

La référence des professionnels